

MICRONIC | SAMPLE STORAGE | CAPS

NOVITAIN CAPCLUSTER-50

DO YOU HAVE SAMPLE COLLECTION TUBES (INCLUDING BLOOD OR URINE SAMPLES) THAT NEED TO BE SECURELY STORED FOR SEVERAL DAYS IN CASE ADDITIONAL ANALYTICAL TESTING NEEDS TO BE DONE? IF SO, THEN MICRONIC OFFERS A UNIQUE POST-ANALYTICAL CAPPING SOLUTION. THE NOVITAIN CAPCLUSTER-50 OFFERS AN IDEAL SOLUTION FOR THE SECURE SHORT-TERM SEALING OF SAMPLE COLLECTION TUBES IN A 50-TUBE FORMAT.

Blood or urine analysis are two of the most important diagnostic tools available to clinicians testing for evidence of infection, disease or abnormalities. Typically blood or urine is collected in vacuum tubes before being sent to a hospital laboratory for testing. After these analytical tests, the vacuum tubes need to be stored for several days in case additional tests need to be performed.

Designed in conjunction with healthcare professionals, Novitain TPE Push Caps offer an ideal solution for the secure short-term sealing of 13.0mm diameter sample collection tubes in a 50-tube format.

These TPE Push Caps are manufactured from a chemically resistant medical grade hydrophobic polymer to ensure stored sample integrity.

Drawing upon proprietary Micronic technology, Novitain Capcluster-50 offers a rapid and economical solution for capping 50 individual sample tubes using the Novitain Manual Capper CV750.

The Novitain TPE Push Caps are available in a Capcluster-50 format. After capping 50 tubes in one turn, the retaining foil can be removed to leave each tube individually capped.

For more information:
e-mail to sales@micronic.com

Micronic Europe B.V.
(Asia, Africa, Europe and Oceania)
Phone: +31 (0)320 277 070
Micronic America LLC
(North, Central and South America)
Phone: +1 484 480 3372

Please visit our website: www.Micronic.com



INFORMATION


- Economical solution for the temporary storage of sample collection tubes
- Caps are provided with a small tab for easy removal
- Caps can be pierced repeatedly without losing seal integrity
- Available in Natural
- Natural TPE caps are made out of Medical approved material (USP Class VI certificate)
- Compatible with sample collection tubes with a diameter of 13.0mm
- Resists chemical solvents




MICRONIC | SAMPLE STORAGE | CAPS

SPECIFICATIONS

MATERIAL	THERMO PLASTIC ELASTOMER
FORMAT	50-TUBE FORMAT (5x10)
RESISTANCE	RESISTS CHEMICAL SOLVENTS
COLORS	NATURAL
AUTOCLAVABLE	YES, INDIVIDUAL CAPS (WITHOUT THE RETAINING FOIL) UP TO TEMPERATURES OF +121°C

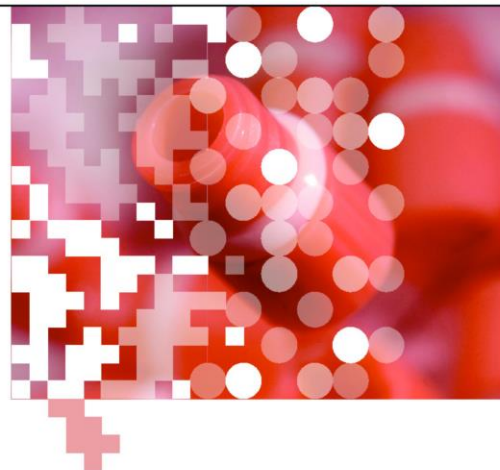

TECHNICAL DRAWING

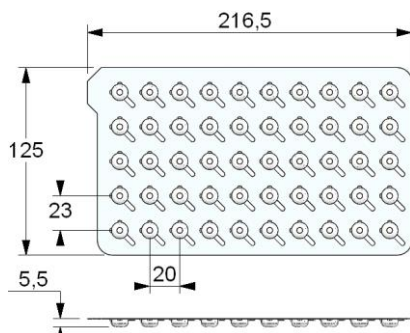

ORDERING INFORMATION

■ Novitain TPE Capcluster for capping 50 individual tubes		
MP35110	Novitain TPE Capcluster Natural for capping 50 individual tubes	box/10

Micronic has a broader product range than the products listed under ordering information. For our complete product portfolio, please send your e-mail to sales@micronic.com or contact us by phone.